

NEVER STOP TO
IMPROVE THE QUALITY AND SERVICE



HOSPITAL CHAIR

VINYL & FABRIC SEATING FOR HEALTHCARE ENVIRONMENTS



Well-experienced



Quality first





Unlimited innovation



Hearted service



Customized service

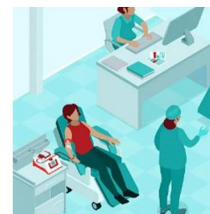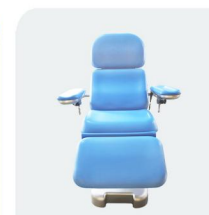

江苏如翼医疗科技有限公司
JIANGSU ROOE MEDICAL TECHNOLOGY CO.,LTD.
Add: Room1311,1313, Guotai New World Plaza, No19,
East Renmin Road,Zhangjiagang, Jiangsu, China.
Tel:+86 512 58810663
Whatsapp/Mobile:+ 86-18962287398;15962399164
Email:info@wanrooemed.com
www.wanrooemed.com www.rooemed.com



江苏如翼医疗科技有限公司
JIANGSU ROOE MEDICAL TECHNOLOGY CO.,LTD.



HOSPITAL CHAIRS

BXD Series Blood Collection Chairs

The BXD series blood collection chairs are specially tailored for medical institutions such as and hospitals blood stations, meeting various medical needs including blood collection, dialysis, blood transfusion, chemotherapy, etc.

BHC001 Series Caregiver Chairs

The BHC001 series caregiver chairs meet the care needs of wards of all levels and provide caregivers with sufficient comfort and care.

BHC003 Series Infusion Chairs

The BHC003 series infusion chairs are exquisitely crafted and offer a comfortable texture, creating an ideal treatment environment for patients.

BHC301 Series Recliners

The BHC301 series recliners are specifically designed for the elderly and patients recovering from surgery. They feature a comfortable texture, adjustable postures, and are easy to move.

BHC302 Elderly Standing Chairs

The BHC302 elderly standing chair not only ensures comfortable sitting and lying positions but also can assist the elderly in standing up through electric control.

CONTENTS

目录

Blood Collection Chairs ————— 01-18

Recliners ————— 19-20

Caregiver Chairs ————— 21-26

Transfusion-Chairs ————— 26-30

Electric Two function Blood Donor Chair

BXD100B PLAY

The Blood donation chair offers better benefits to individuals such as volunteers who donate blood because it can also be used as a dialysis chair, a medical and rehabilitation chair, and even a bed. Electric Blood donation chair offers adjustable and flexible sitting surfaces in addition to its already impressive technical features and specifications. Some of the features of Blood Chair include a high quality steel frame, multiple functional options, and maximum load bearing capacity totaling 200Kg.



Seat position



Semi-fowler position



Horizontal position



CPR trendelenburg position



Manual CPR



Dining table

Product model	BXD100B
Chair open size	1900*850*560mm
Chair body size	920*850*560mm
Package size	1510*800*610mm/72kg
Seat width	600mm
Arm rest adjustment height	none
Height Adjustable	none
Back adjustment Angle range	-10° ~ 75°±5°
Leg plate Adjustment Angle range	-80° ~ 12°±5°
Trendelenburg position	-12°
IV pole from clamp to top	860~1500mm
IV pole from floor to top	1260~1900mm
Frame Material	Steel
Matress Material	PVC
Wheel	Four 5-inch wheels with brakes 4-inch optional
Max loading	200kg
N.W/G.W.	62kg/75kgs
Power Supply	AC220V 50Hz/60Hz

1. Designed according to the principle of human body engineering, luxurious, comfortable and durable. You won't get tired if you sit too long.
2. The blood collection chair consists of four main parts: backrest, seat, leg plate and foot pedal.
3. It carries out the linkage movement of back, legs and back legs through the electric remote control.
4. Wide and comfortable handrails support blood collection.
5. Extendable foot pedal boards provide patient comfort and comfort.
6. Q235 steel pipe is durable, cushion is made of imported PVC high-quality artificial leather, polyurethane one-time molding.
7. Four 5-inch casters, each wheel can be individually locked and the seat can be fixed in the desired position.
8. The chair can be easily used for blood collection and dialysis procedures.
9. Optional: Electric CPR, IV Stand, Oxygen Tank Stand, Central locking castors, paper holder.

Manual Reclining Blood Draw Chair

BXD100B Basic PLAY

Our Reclining Phlebotomy chair can be adjusted by a handle next to the seat for examination of the patient. The blood collection chair uses a gas spring mechanism instead of electric adjustment. Adjustable Phlebotomy Chair's wheels move the chair from one place to another. The maximum safe bearing capacity of 200kg also provides high stability in the horizontal position. Phlebotomy draw chairs are widely used in blood donation centers and medical centers.



Seat position



Semi-fowler position



Backrest adjustment, backleg linkage



Horizontal position



Cpr trendelenburg position



Manual adjustable body position

Product model	BXD100B Basic
Chair open size	1900*850*560mm
Chair body size	920*850*560mm
Package size	1620*820*600mm/72kg
Seat width	600mm
Arm rest adjustment height	none
Height Adjustable	none
Back adjustment Angle range	-10° ~ 75°±5°
Leg plate Adjustment Angle range	-80° ~ 12°±5°
Trendelenburg position	-12°
IV pole(from clamp to top)	860~1500mm
IV pole(from floor to top)	1260~1900mm
Paper holder and oxygen holder diameter	none
Frame Material	Steel
Matress Material	PVC
Wheel	Four 5-inch wheels with brakes 4-inch optional
Max loading	200kg
N.W/G.W.	62kg/75kgs

1. Low carbon steel solid structure, appearance of high strength fiber material and multi-layer baking paint, with artistic sense, modern sense, durability and easy to clean.
2. The mattress is made of high-quality medical polyurethane foam forming, which is easy to clean and maintain. Its antistatic and antifouling performance is superior to other similar products.
3. The combined backrest and legrest movement promotes a relaxed reclining/lying flat position through gas spring assisted movement.
4. Release the handle to quickly adjust the seat position. The release handle can be operated from either side of the recliner.
8. Super wide upholstered handrails maximize comfort. Optional dining table board.
9. Manual CPR is standard.
10. It is an ideal blood collection choice for blood stations, centers and blood donation rooms.
11. Human-oriented original products can give blood donors a sense of security and comfort, relieve psychological pressure, and make the whole process relaxed and happy.
9. Optional: IV Stand, Oxygen Tank Stand, Central locking castors, paper holder.

Electric Three function Blood Donor Chair

BXD100B Plus PLAY

The adjustable phlebotomy chair has an adjustable seat height and is designed to meet the various healthcare needs of patients for examination, diagnosis and treatment. Electric phlebotomy chairs have a high-density foam cushion that is flexible enough not to cause postural strain when sitting or lying for long periods of time, and PVC leather that is durable and luxuriously comfortable. The dialysis recliner chairs' retractable, foldable footrests provide comfort for patients during dialysis.



1. Steel frame coated with epoxy powder with high mechanical strength, chemical and corrosion, beautiful, modern, durable and easy to clean.
2. The mattress is made of high-quality medical polyurethane foam forming, which is easy to clean and maintain. Its antistatic and antifouling performance is superior to other similar products.
3. With its easy-to-use controller, the seat can stop in any intermediate position and the patient can be positioned comfortably.
4. The hand controller adopts easy-to-read button design, can automatically adjust the sitting position or flat bed, and can individually adjust the seat Angle, seat height, backrest Angle, pedal Angle, CPR emergency headdown positioning (or trendelenburg position).
5. Four 5-inch silent medical casters with independent brakes.
6. Removable headrest The pillow can be adjusted according to the height of each patient.
7. Cushioned armrest for comfortable sitting position.
8. Electric power control and adjustment of sitting position and foot position to create a safe and comfortable environment for blood donors.
9. There are holes for brine rack around, which is convenient to plug and install the brine rack in the plug pipe.
10. Ideal for dialysis, blood collection, respiratory and cardiac care.
11. Optional: Electric CPR, IV Stand, Oxygen Tank Stand, Central locking castors, paper holder.



Seat position




Horizontal position



CPR Trendelenburg position



Double-layer rotatable armrest

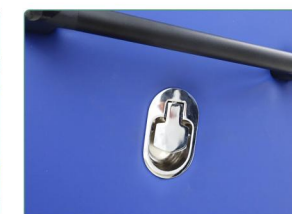

Manual CPR



Dining table

Product model	BXD100B Plus
Chair open size	1900*850*560/820mm
Chair body size	920*850*560/820mm
Package size	1610*820*600mm/72kg
Seat width	600mm
Arm rest adjustment height	none
Height Adjustable	560~820mm
Back adjustment Angle range	-10° ~ 75°±5°
Leg plate Adjustment Angle range	-80° ~ 12°±5°
Trendelenburg position	-12°
IV pole (from clamp to top)	860~1500mm
IV pole (from floor to top)	1260~1900mm
Frame Material	Steel
Mattress Material	PVC
Wheel	Four 5-inch wheels with brakes □4-inch optional□
Max loading	200kg
N.W/G.W.	62kg/75kgs
Power Supply	AC220V 50Hz/60Hz

Electric Three function Blood Donor Chair

BXD100B Pro PLAY

The nursing recliner makes patient care, treatment and recovery more comfortable and provides greater convenience for caregivers in a variety of clinical Settings. Helps your patients relax while improving efficiency.



1. Fully electric and accurately adjustable blood drawing chair, which can easily locate blood donors.
2. Electric remote control adjustment: overall lifting, back lift, leg lift, back leg linkage.
3. Double armrest design, the upper armrest can be adjusted left and right, to avoid the blood donor arm sliding and unstable placement.
4. Velcro headrest for head and neck comfort during blood donation.
5. 4 Casters with brakes for easy movement.
6. Foldable foot pedals to provide adequate foot support and prevent the patient from sliding down.
7. Manual CPR to assist in emergency situations.
8. Optional: Electric CPR, IV Stand, Oxygen Tank Stand, Central locking castors, paper holder.



Seat position




Backrest adjustment, backleg linkage



Horizontal position



CPR Trendelenburg position

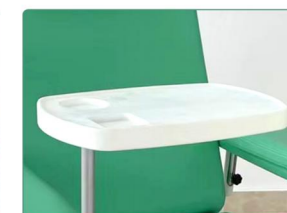

Dining table



Movable upper handrail

The reclining phlebotomy chair can be powered upright, partially reclined or fully reclined. With synchronized headrest and leg support for easy adjustment and quick recovery of positioning. With an upholstered seat and a back with adjustable armrests, patients are sure to stay comfortable. The electric phlebotomy chair keeps patients comfortable and can better accommodate patients who may faint.

Product model	BXD100B Pro
Total Length	1950mm+20mm
Seat width	600mm (total width 900mm including armrest)+20mm
Seat height	580~780mm+20mm
Backrest length	860mm+20mm
Height of armrest	230~290mm+20mm
Castor	4xφ 100mm swivel castors with separate brakes
Pillow size	410mm*230mm*80mm+20mm
Backrest Adjustment	(-12°~75°)+5°
Largest Adjustment	(-75°~12°)+5°
Weight	70KGS+3KGS
Safe Maximum Load	240KGS

Head Adjustable Phlebotomy Chair

BXD101 PLAY



- 1.Usage: Ideal for blood collection, dialysis, chemotherapy in medical units.
- 2.Structure: 4-section (headrest, backrest, seat, leg rest) electric blood collection chair. The seat maintains a comfortable angle of about 5 degrees. The Plus version has 5 sections with an additional footrest.
- 3.Frame: The base, backrest cover and supports are made of ABS plastic by molding, which is aesthetic, dustproof, durable and easy to clean.
- 4.Control: Each part can be adjusted in height and angle via electric buttons to meet different patients' and postures' needs.
- 5.Surface: Covered with seamless PU leather, fireproof and anti-static. Thick sponge conforms to the body shape.
- 6.Armrest: Adjustable in height, position and angle. Anti-slip arc ensures arm stability.
- 7.Configuration: Equipped with a push rod to assist patient transfer. Casters with brakes for easy fixing.
- 8.Optional: Multimedia bracket to add treatment fun; footrest and posture memory; foot control for quick posture adjustment.
- 9.Optional: Seat memory function to reach the preset position with one click.


BXD101PLUS



Product model	BXD101	BXD101PLUS
Chair dimension	1710±60mm×560±40mm	1710±60mm×560±40mm
Height	600±40mm—830±60mm	600±40mm—830±60mm
Back section adjust range	0°~75°±5°	0°~75°±5°
Seat section angle	5°±2°	5°±2°
Leg section adjust range	10°~65°±5°	10°~65°±5°
Arm rest tilt angle	-15°±5°~15°±5°	-15°±5°~15°±5°
Arm rest swing out	100±20mm	100±20mm
Arm rest height adjustment	0-12°	0-12°
Arm board slide	100±20mm	100±20mm
Foot section adjust range	none	85°± 5°
Foot section telescope	none	100mm± 20mm
Power input	220V±10% 50Hz	220V±10% 50Hz
Safety working load	200kg	200kg
Multimedia support arm	optional	optional

Electric Blood Collection Chair

BXD130 PLAY

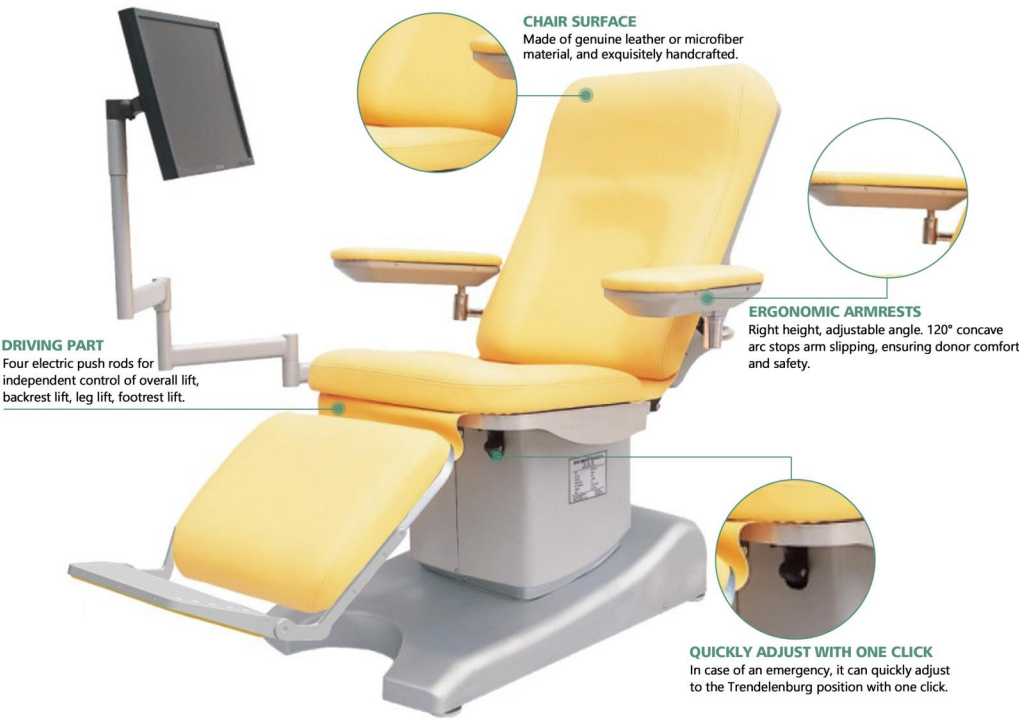

- 1.Usage: Ideal for blood collection, dialysis, chemotherapy in medical units.
- 2.Structure: 4-section (headrest, backrest, seat, leg rest) electric blood collection chair. The seat maintains a comfortable angle of about 5 degrees. The Plus version has 5 sections with an additional footrest.
- 3.Frame: The base, backrest cover and supports are made of ABS plastic by molding, which is aesthetic, dustproof, durable and easy to clean.
- 4.Control: Each part can be adjusted in height and angle via electric buttons to meet different patients' and postures' needs.
- 5.Surface: Covered with seamless PU leather, fireproof and anti-static. Thick sponge conforms to the body shape.
- 6.Armrest: Adjustable in height, position and angle. Anti-slip arc ensures arm stability.
- 7.Configuration: Equipped with a push rod to assist patient transfer. Casters with brakes for easy fixing.
- 8.Optional: Multimedia bracket to add treatment fun; footrest and posture memory; foot control for quick posture adjustment.
- 9.Optional: Seat memory function to reach the preset position with one click.

BXD170 PLAY



Product model	BXD130	BXD170
Dimension (L*W*H) (mm)	1950*580*550-700	2000*580*550-750
Back section adjustable	8°-70°	8°-70°
Leg section adjustable	10°-60°	20°-75°
Trendelenburg/reverse trendelenburg	8°-15°	8°-15°
Dimension (mm)	500*800	480*160
Swing out	0°-45°	0°-45°
Tilt angle	-10°	100mm
Slide distance (mm)	100mm	150mm
Table top size (mm)	380*280	450*320
Power supply	AC220V/50HZ	AC200V/50HZ
Working voltage	24V (electric type)	24V (electric type)
Capacity	more than 220kg	more than 240kg

Electric Dialysis Bed

BXD188B PLAY



1. Product Use: It is suitable for the collection of component blood in blood supply stations, as well as blood collection, dialysis, infusion chemotherapy, etc. in hospitals, fully meeting the needs of medical institutions.
2. Product Structure: It has a four-section structure (backboard, seat board, leg board, footrest board) with electric adjustment, allowing for easy changes in body position to ensure convenient operation and patient comfort.
3. Main Frame: The bed frame is made of high-quality steel, treated with electrostatic plastic spraying at high temperature, which is anti-corrosive and wear-resistant. The backboard and base cover are made of ABS engineering plastics, and aluminum alloy lifting guardrails are assembled on the side, ensuring outstanding quality.
4. Control System: It is electrically controlled and driven by electric push rods, enabling precise adjustment of the lifting of the back and legs, and easily achieving the best blood collection or treatment posture.
5. Surface Material: The chair surface is covered with PU leather, which is seamless, fireproof, and anti-static. It is filled with thickened sponge that perfectly conforms to the body shape of each patient, providing the best relaxation effect and improving the treatment experience.
6. Armrest Design: Based on ergonomic principles, the height and angle conform to the natural stretching curve of the human arm, ensuring a comfortable sitting posture and preventing the arm from sliding.
7. Other Functions: It can be converted into a temporary hospital bed when laid flat. ABS universal casters are equipped at the bottom, with a load-bearing capacity of up to 240 kg, ensuring quiet and stable movement.

This medical bed is perfect for blood collection, dialysis, infusions, and chemotherapy in hospitals. It features a four-section electric-adjustable design for comfort and ease of use. The electric actuator ensures precise positioning for treatments. Made with soft microfiber and durable ABS components, it offers patient safety with an aluminum guardrail and non-slip features. The high-quality steel frame is corrosion-resistant, and the bed is equipped with smooth-rolling ABS castors for easy movement in any medical setting.

Product model	BXD188B
Dimension(L*W*H)(mm)	1950*780-900*520
Back section adjustable	0°~ 70°
Leg section adjustable	0°~ 60°
Dimension (mm)	500*130
Swing out	0°~ 45°
Slide distance (mm)	0-100
Power supply	AC220V/50HZ
Working voltage	24V (electric type)
Capacity	more than 220kg

Blood Drawing Chair

BXD103 PLAY

- 1. Welded powder coated frame is used to increase stability and for better support.
- 2. The seat has a PU cushion which allows doctors to draw blood or measure blood pressure more easily.
- 3. Height and depth adjustable flip arm can be locked in the right position.
- 4. The arm allows easier access to the patient when the doctor has to draw blood or measure blood pressure.
- 5. Optional: Drawers - Drawers can be placed on either side of the chair for easy access to phlebotomy tools and equipment.



Product model	BXD103
Overall size(L*W*H)(mm)	680*650*950
Back section size(mm)	540*390
Seat size(mm)	550*520
Side armrest size(mm)	387*187*40
Front armrest size(mm)	465*160
Capacity	more than 180kg

Function Blood Draw Chair

BXD104

- 1.Product Usage: Ideal for blood collection, dialysis, chemo chairs in medical units.
- 2.Product Structure: Four - section (back, seat, leg, footrest) electric blood collection chair, adaptable to various postures.
- 3.Main Frame: Steel frame for high load - bearing. ABS - molded backrest cover, nice - looking and easy to clean.
- 4.Control System: Backrest and legs adjusted via manual controller. Adjustable armrests (height, position, angle), electric operation for multiple postures.
- 5.Surface Material: PU - covered seat, seamless, fire - resistant, anti - static. Thick sponge for comfort.
- 6.Armrest Design: Concave armrests prevent slipping. Stay horizontal during backrest adjustment.
- 7.Other Functions: Equipped with push rod and caster wheels with brakes for stability and mobility.



Product model	BXD104
Dimension(L*W*H)(mm)	1915±60*560±40*540±40
Head section adjustable	45°±5°
Back section adjustable	30°±5°
Seat section adjustable	11°±2°
Leg section adjustable	85°±5°
Swing out	45°±5°
Tilt angle	-15°±5°~15°±5°
Slide distance (mm)	100±20mm
Power supply	AC220V/50HZ
Capacity	more than 135kg

Manual Blood Donor Chair

BXD105 PLAY

- 1.Usage: It's an ideal medical chair for blood collection, dialysis, chemotherapy in medical units.
- 2.Structure: This four - section blood collection chair (backrest, seat, leg board, optional footrest) has a main frame of high - quality steel with powder coating.
- 3.Control: The chair has a control crank for manual adjustment of back and leg angles, keeping patients comfortable during treatment.
- 4.Surface: PU leather, seamless, fireproof, anti - static. Thick sponge for body fit.
- 5.Armrest: Adjustable for proper arm positioning.
- 6.Options: Mobile wheels, saline stands, oxygen brackets, paper holders, footrests available.



Product model	BXD105
Dimension(mm)	1650(curve length)* 560mm
Back section adjustable	0-70°
Leg section adjustable	0-60°(leg section and back section work as one unit)
Armrest Dimension(mm)	450*160
Armrest Height adjustment(mm)	100mm
Armrest outward swing angle	0-45°
Arm board concave radian	R108
Package (1pc/wooden case)(mm)	1400*750*440(mm)
N.W.	50.5kg
G.W.	66kg

BXD105 Basic PLAY

1. Use: Ideal medical chair for blood collection, dialysis, chemo in med units. Made of high - strength fiber, multi - layer paint.
2. Structure: Three - section (backrest, seat, leg board) integrated blood collection chair, ergonomic design for comfortable sitting.
3. Main Frame: Made of high - quality steel with powder - coated rigid frame, strong stability and high load - bearing.
4. Surface Material: Chair surface covered with seamless, fireproof and anti - static PU leather, filled with high - density sponge for a comfortable touch.
5. Control System: Manual adjustment of backrest and leg rest with gas spring assistance for easy tilt/reclining.
6. Armrest Design: Swing - adjustable to different angles and height - adjustable to match patients' postures.



Product model	BXD105 Basic
Dimension(L*W*H)(mm)	1650(curve length)* 560
Back section adjustable	0-70°
Leg section adjustable	0-60 °(leg section and back section work as one unit)
Armrest Dimension(mm)	480*160
Armrest Height adjustment(mm)	54-82
Armrest outward swing angle	0-45 °

Elderly Chair

BHC301A PLAY

1. Covered with PVC calendered artificial leather, filled with sponge, widened armrests, expandable leg rests, and comfortable feel.
2. The overall frame is constructed with high-quality steel sprayed plastic pipes, which has strong stability.
3. The backrest is equipped with a handle, allowing the chair to be easily moved.
4. It is equipped with a lockable small table board. There is a circular groove on the small table board for placing cups.
5. With five-inch casters, it moves smoothly and quietly.
6. The mechanical structure is manually operated and can be adjusted to three different angles.



BHC301B PLAY

1. It is covered with PVC calendered artificial leather and filled with shaped sponge. The armrests are widened, offering a comfortable touch.
2. The whole frame is built with high-quality steel pipes coated with sprayed plastic, which ensures strong stability.
3. The backrest is equipped with a handle, enabling the chair to be moved easily.
4. It is equipped with a lockable small table board, on which there is a circular groove for placing cups.
5. With five-inch casters, it moves smoothly and quietly.
6. The mechanical structure is manually operated and can be adjusted to three different angles.
7. A footrest is provided to make the sitting posture more comfortable.



BHC301C PLAY

1. Covered with PVC calendered artificial leather, filled with sponge, widened armrests, expandable leg rests, and comfortable feel.
2. The entire frame is constructed with high-quality steel pipes coated with sprayed plastic, ensuring strong stability.
3. There are two small table boards attached to the armrests on both sides, making it more convenient to place items.
4. Equipped with five-inch casters, it moves smoothly and quietly.
5. It is equipped with a gas spring to control the angle of the backrest, which can be adjusted to three different angles.
6. Can be quickly adjusted to the shock position to deal with emergencies.



Product model	BHC301A	BHC301B	BHC301C
Product size	1100*660*1180mm	1400*730*1190mm	1080*770*1280mm
Packing Size	950*700*720mm	960*800*740mm	960*800*740mm
N.W.(kg)	41	45	44
Cushion size	530*510mm	530*510mm	530*510mm
Cushion height	520mm	520mm	520mm
Plate size	650*380mm	650*380mm	/

Electric standing-up chair

BHC302 PLAY

The elderly recliner provides a more comfortable treatment experience for the elderly and patients. The surface of the electric elderly chair is covered with genuine leather and filled with thickened sponge, ensuring a comfortable sitting posture and providing sufficient support for the spine and neck. It is equipped with a remote control, allowing patients to independently operate and change the angle and posture of the chair, as well as adjust their resting position. Meanwhile, it can tilt forward as a whole to assist those with limited strength in standing up, and it is equipped with casters for convenient movement. (The functions of waist warmth and massage are optional.)



Product model	BHC302
Overall size(L*W*H)	1030*750*1040mm
Seat size(W*D*H)	550*580*460mm
Full recliner length	1840mm
Back size	60*70cm
Lifting height	57-70cm
Back section adjustment	105°-171°
Foot section adjustment	9°-20°
Lifted angle	0 -20°
Loading capacity	100kg

Luxurious Noon Break Chair

BHC001C PLAY

1. Aluminum alloy armrests with black crusty soft padding.
2. High - quality bright carbon steel tubes, electrostatically sprayed with imported powder and baked at high - temp. Corrosion - resistant, coating won't peel.
3. Seat cushion surface in PVC calendered artificial leather for great tactile feel.
4. Thick - padded chair, with head area thickened when unfolded for comfy sit - lie.
5. Silent casters for quiet, smooth, jam - free movement.
6. Telescopic springs inside for free extension and retraction, easily convertible to a single bed.



BHC001N

1. It features a sofa - style appearance design, with hidden casters and a steel frame structure, presenting an elegant and grand look.
2. The surface of the seat cushion is made of PVC calendered artificial leather, offering an outstanding tactile experience.
3. Telescopic springs are equipped inside. During the unfolding process, it can self - position, extend and retract freely, and can be easily transformed into a single bed.



Product model	BHC001C	BHC001N
Product size	780*650*450mm	830*850*420mm
Expanded size	2060*650*450mm	2100*850*620mm
Backrest height	1060mm	960mm
Packing size	1280*660*350mm	850*850*700mm
N.W.(kg)	38	55

BHC001O & BHC001O(Widened) PLAY

1. It adopts a sofa - style appearance design, combined with hidden casters and a steel frame structure, presenting an elegant and generous look.
2. The entire surface is covered with PVC calendered artificial leather, providing an outstanding tactile experience.
3. There is a hidden storage cabinet under the seat cushion, which has a strong storage function.
4. Telescopic springs are equipped inside. During the unfolding and folding process, it can self - position, be extended and retracted freely, and is easy and convenient to operate.
(A widened version is available.)



(A widened version is available.)

Product model	BHC001O	BHC001O(Widened)
Product size	880*740*420mm	850*860*450mm
Expanded size	2000*740*640mm	2000*860*640mm
Backrest height	900mm	900mm
Packing size	900*760*920mm	870*880*920mm
N.W.(kg)	60	60

Luxurious Accompanier's Chair

BHC001F PLAY

- 1. The surfaces of the seat cushion and the backrest are covered with PVC material through calendaring imitation leather, offering an excellent tactile experience.
- 2. The steel frame is welded with spray-painted steel pipes, and the armrests are covered with sponge.
- 3. There are telescopic springs inside, enabling easy retraction and extension.
- 4. After unfolding, the cushion at the head part is thickened, providing a comfortable lying experience.



BHC001K PLAY

- 1. The overall surface is covered with PVC material through calendaring imitation leather, providing an excellent tactile experience.
- 2. The concealed steel frame structure is made of spray-painted steel pipes, which are aesthetically pleasing and grand.
- 3. There are telescopic springs inside, allowing for easy retraction and extension.
- 4. The head part has an increased height design after unfolding, ensuring a comfortable lying experience.



Product model	BHC001F	BHC001K
Product size	660*680*410mm	800*710*430mm
Expanded size	1860*680*400mm	2040*710*430mm
Backrest height	930mm	980mm
Packing size	700*700*130mm/ 780*560*280mm	780*60*350mm/ 740*740*180mm
N.W.(kg)	32	36

BHC001M

- 1. It has a luxurious sofa-style appearance, equipped with hidden casters and a steel frame structure.
- 2. It is aesthetically pleasing and grand, and the solid wood armrests offer an excellent texture.
- 3. It is equipped with telescopic springs inside. During the unfolding process, it can self-position, and it can be easily retracted and extended.
- 4. It can be easily transformed into a single bed.



BHC001L

- 1. It features a specially widened design, and there is also a backrest after it is unfolded.
- 2. The fabric is linen, which takes into account both comfort and cost-effectiveness.
- 3. There are telescopic springs inside, allowing for easy retraction and extension.



Product model	BHC001M	BHC001L
Product size	880*790*420mm	800*1020*430mm
Expanded size	2000*790*610mm	2080*1020*430mm
Backrest height	990mm	980mm
Packing Size	900*810*1010mm	920*800*430mm
N.W.(kg)	52	55

Luxurious Accompanier's Chair

BHC001 PLAY

- 1. The steel frame is welded with spray-painted steel pipes, and the fabric is made of PVC calendered imitation leather.
- 2. The armrests are covered with plastic parts to improve the tactile sensation, and other materials can also be selected for covering.
- 3. A locking device can be optionally equipped to lock the deformed state of the chair to ensure safety.
- 4. A gas spring can also be optionally equipped to make the retraction and extension of the chair more smooth and effortless.



BHC001A

- 1. It has a stowable design, which is convenient for storage and saves space.
- 2. The steel frame is welded with spray-painted steel pipes, and the fabric is made of PVC calendered imitation leather.



Product model	BHC001	BHC001A
Product size	605*630*420mm	710*640*300mm
Expanded size	1900*630*580mm	1900*640*200mm
Packing size	1190*620*220mm	730*650*320mm
N.W.(kg)	23	21

BHC001S

- 1. The steel frame is welded with stainless steel pipes, which is aesthetically pleasing and grand, and has a long service life.
- 2. The fabric is PVC calendered imitation leather.
- 3. The armrests are covered with solid wood to improve the tactile experience.



BHC001G

- 1. Made of flannel fabric, it is soft, delicate and provides a comfortable tactile experience.
- 2. The iron frame is welded with spray-painted steel pipes.
- 3. After unfolding, the head part can be lifted up, enhancing the comfort.
- 4. It is convenient to retract and extend, and easy to store.



Product model	BHC001S	BHC001G
Product size	650*660*420mm	750*240*980mm
Expanded size	1900*660*580mm	1860*750*330mm
Packing size	1190*620*220mm	700*240*920mm
N.W.(kg)	18	18.5

Luxurious Transfusion-Chair

BHC003A

This meticulously crafted luxury infusion chair is designed to comprehensively meet these demands of patients, starting from every single detail. It integrates advanced design concepts and high - quality material craftsmanship, aiming to bring patients an unprecedented and almost home - like comfortable infusion experience, enabling every user to feel attentive care and solicitude during the infusion process.



- 1.The chair uses top - grade PVC calendered artificial leather, mimicking real leather for durability. Filled with thick, high - density sponge, it eases long - term sitting discomfort. A three - angle adjustment lets users find the best position for a comfy infusion.
- 2.There's a foldable table on the side for meds, drinks, etc. and a sturdy metal storage basket for personal items, keeping the area neat.
- 3.Its frame is made of high - quality spray - coated steel pipes, resistant to rust and looking modern. Widened solid wood armrests offer more support and add elegance.
- 4.An electric adjustment option is available. With a simple button press, it's easy and precise, especially helpful for patients with limited mobility, enhancing comfort and accessibility.

Luxurious Gas Spring Transfusion-Chair

BHC003F

- 1. This product features a cozy sofa - style design. Its top - notch PVC calendered artificial leather surface closely mimics real leather, being both high - end and durable, resistant to scratches and suitable for daily use in various settings.
- 2. Hidden within is a well - engineered steel frame. This concealed design maintains the sofa's sleek appearance while providing robust support. Made of high - quality steel, the frame ensures stability, withstanding substantial weight without deforming and adding elegance to any room.
- 3. Inside the sofa, a state - of - the - art gas spring is ingeniously installed. It enables easy adjustment of the chair to any angle for reading, napping or lounging, offering smooth and precise angle changes for a personalized, more comfortable and functional seating experience.
- 4. At the bottom are four high - quality silent casters with noise - reducing technology. They allow for quiet and easy movement of the sofa in the living space, whether for furniture rearrangement, cleaning or changing the room layout, adding daily convenience.



Product model	BHC003A	BHC003F
Product size	900*750*1100mm	900*850*635mm
Packing size	980*600*440mm	870*850*635mm (Semi-assembled)
N.W.(kg)	31	38
Backrest height	/	1050mm

Transfusion-Chair

BHC003

1. The fabric is made of PVC calendered artificial leather, filled with thickened sponge, providing a comfortable texture.
2. The whole product can be adjusted to three different angles, ensuring a comfortable posture during infusion.
3. It is welded with spray-coated steel pipes. The armrest part is designed with widened solid wood armrests, and a metal storage basket is equipped on the side.



BHC003B

1. The product features a covering crafted from PVC calendered artificial leather, which is filled with thickened sponge. This combination imparts an exceptionally comfortable texture.
2. Its structure is constructed through the welding of spray-coated steel pipes, ensuring durability.
3. At the armrest area, widened solid wood armrests have been ingeniously designed, enhancing both comfort and aesthetics.
4. On the side, a metal storage basket is conveniently equipped, adding to the product's functionality.

Product model	BHC003	BHC003S	BHC003B	BHC003D
Product size	770*680*1080mm	770*680*1080mm	770*680*1080mm	820*700*1020mm
Packing size	920*620*355mm	920*620*355mm	820*620*310mm	900*560*310mm
N.W.(kg)	21.5	22	19	25

BHC003S

1. The frame, welded with stainless steel pipes, resists dirt and corrosion.
2. The fabric is made of PVC calendered artificial leather, filled with thickened sponge, offering a comfortable texture.
3. The whole can be adjusted to three different angles, ensuring a comfortable posture during the infusion process.
4. A stainless steel storage basket is equipped on the side, and a foldable small table board on the side can be optionally equipped.



BHC003D

1. The fabric is made of PVC calendered artificial leather, and it is welded with spray-coated steel pipes.
2. Widened solid wood armrests are designed at the armrest part, and a metal storage basket is equipped under the seat.
3. A children's version is available, and a small stool is provided for convenient operation.

